


patient commission/configuration ID _____



Item	12mm	Description	Potential Problem	Measure	Recommended Inspection, Potential Replacement*	Latest Replacement after
System Joint						
1	SB7049-L0590	bearing nut	wear	replacing	every 6 months	36 months
3	SF0391-01	feather key with pin	breakage	replacing	every 6 months	36 months
5	GS1407-040	sliding washer, Ø = 14mm	wear	replacing	every 6 months	18 months
5	GS1407-045					
5	GS1407-050					
5	GS1407-055					
5	GS1407-060					
6	SH0461-L/AL	cover plate, left lateral or right medial**	wear	replacing	every 6 months	36 months
6	SH0461-R/AL	cover plate, left medial or right lateral**				
7a	VE3771-050/10	O-ring for securing the alignment screw**	wear	replacing	every 6 months	18 months
7c	FE1611-01	coil spring**	wear	replacing	every 6 months	18 months
7d	GS0705-525	sliding bushing**	wear	replacing	every 6 months	18 months
7e	PN0004-L02	stop damper**	wear	replacing	every 6 months	18 months
7f	SH0491-01	plunger**	wear	replacing	every 6 months	36 months
8a	VE3771-08/15	O-ring for securing the alignment screw**	wear	replacing	every 6 months	18 months
8c	VE3771-04/10	O-ring for securing the motion limiting screw**	wear	replacing	every 6 months	18 months
9	SC1404-L10	countersunk flat head screw with hexalobular socket**	wear	replacing	every 6 months	36 months
10	SC1404-L10	countersunk flat head screw with hexalobular socket (axle screw)**				
6-10	SH4971-L/AL	functional unit plug + go modularity, left lateral or right medial***	loss of function, wear	replacing	every 6 months	36 months
6-10	SH4971-R/AL	functional unit plug + go modularity, left medial or right lateral***				

patient commission/configuration ID _____

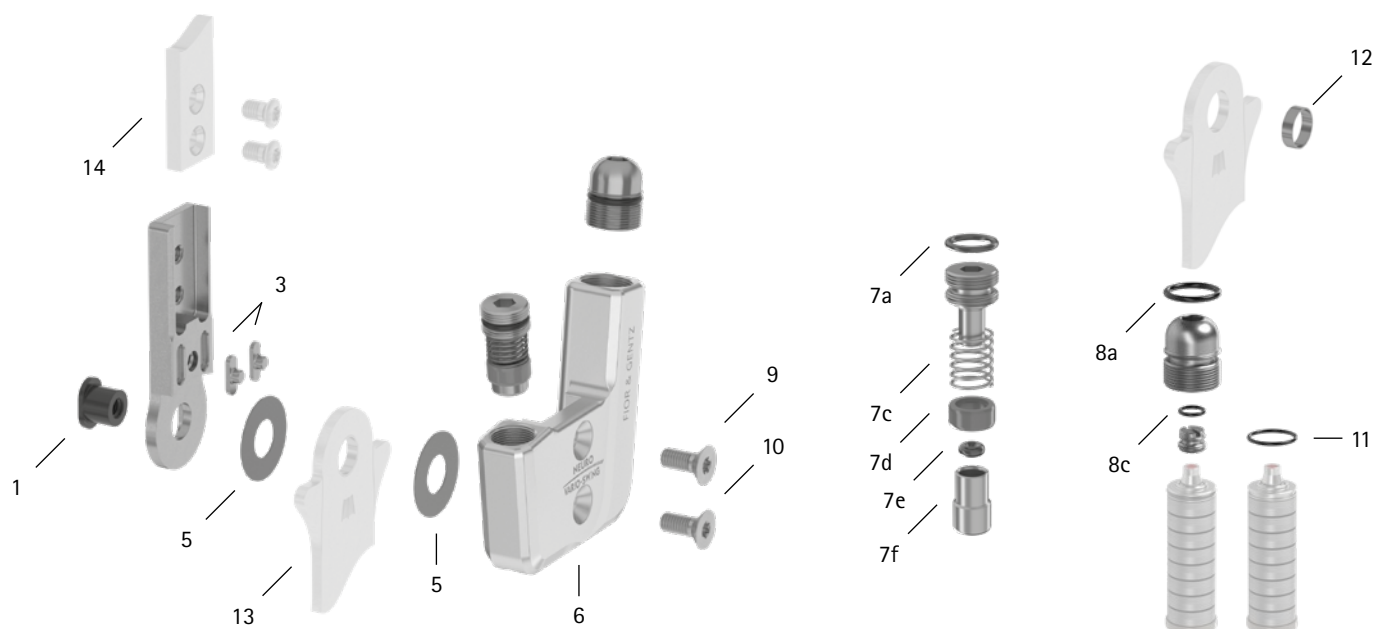

Item	12mm	Description	Potential Problem	Measure	Recommended Inspection, Potential Replacement*	Latest Replacement after
Spring Units						
11	VE3771-06/10	O-ring for securing the spring unit	wear	replacing	every 6 months	18 months
-	SF5801-15/03	spring unit, blue, normal, max. 15° range of motion	wear	replacing	every 6 months	18 months
-	SF5801-15/06	spring unit, green, medium, max. 15° range of motion		realigning disc springs with pliers		
-	SF5801-10/12	spring unit, white, strong, max. 10° range of motion	radial move of disc springs	greasing spring unit with spray oil (article no. FT3000-15)		
-	SF5801-10/19	spring unit, yellow, very strong, max. 10° range of motion	noise of spring unit			
-	SF5801-05/33	spring unit, red, extra strong, max. 5° range of motion				
System Stirrups, System Anchors and System Side Bars						
12	BR0807-L020	sliding bushing	wear	replacing	every 6 months	18 months
13	see catalogue	system stirrup	wear or breakage	replacing	every 6 months	48 months
14	see catalogue	system anchor and system side bar	wear or breakage	replacing	every 6 months	if required

* depending on the assessment of the distributor of the custom-made product regarding the patient's usage behaviour
** is part of the functional unit
*** included system components can be exchanged separately