


patient commission/configuration ID _____

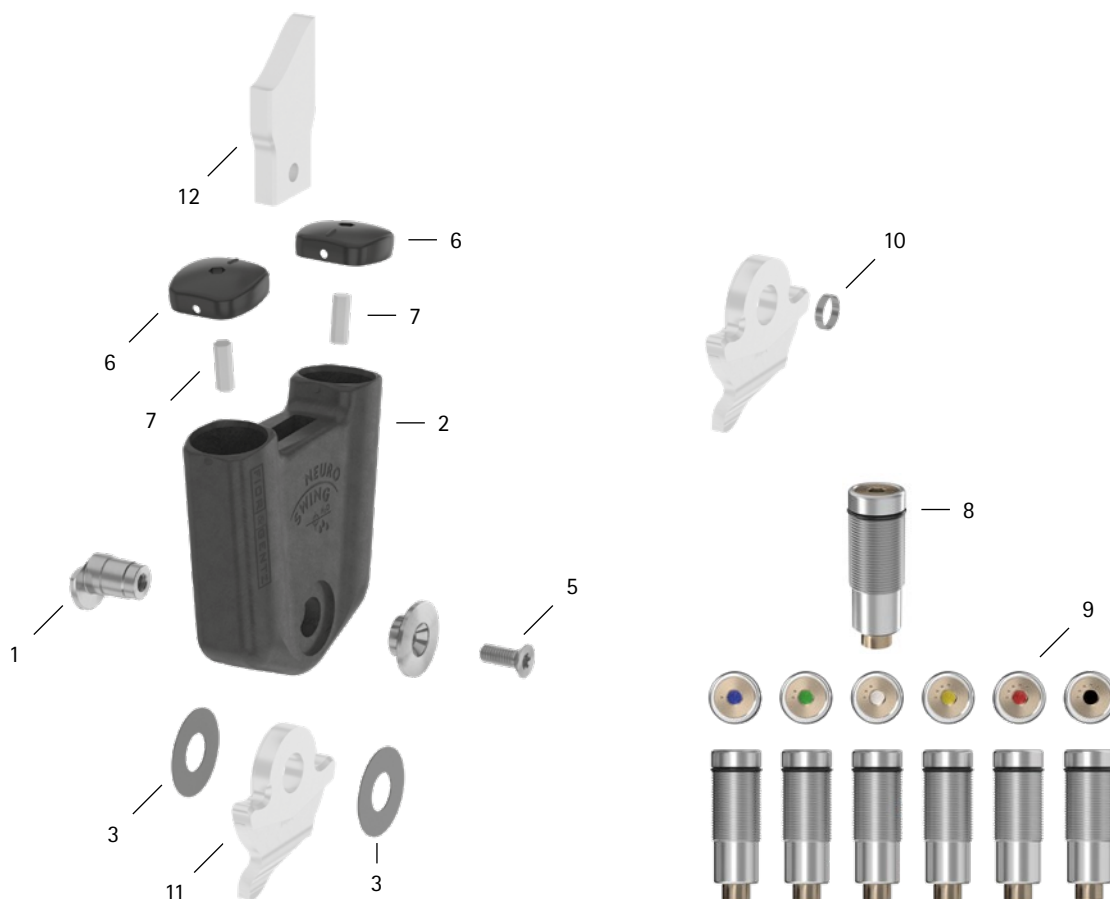

Item	24mm	Description	Potential Problem	Measure	Recommended Inspection, Potential Replacement*	Latest Replacement after
System Joint						
1	SF0597-C/1	bearing nut	wear or breakage	replacing	every 6 months	48 months
2	SF0507-C	joint case	wear or breakage	replacing	every 6 months	48 months
3	GS3013-040	sliding washer, Ø = 30mm	wear	replacing	every 6 months	18 months
3	GS3013-045					
3	GS3013-050					
3	GS3013-055					
3	GS3013-060					
5	SC1416-L16	countersunk flat head screw with hexalobular socket	wear	replacing	every 6 months	36 months
6	SF0597-C/4/1	spring unit cover	wear	replacing	every 6 months	if required
7	SF0597-C/4/0	hexagonal pin	wear	replacing	every 6 months	if required
Spring Units and Static Dorsiflexion Stop						
8	VE3771-180/15	O-ring for securing the spring unit and static dorsiflexion stop	wear	replacing	every 6 months	18 months
9	SF5807-C/15/12	spring unit, blue, normal, max. 15° range of motion	noise of spring unit/ static dorsiflexion stop	replacing	every 6 months	18 months
9	SF5807-C/15/25	spring unit, green, medium, max. 15° range of motion				
9	SF5807-C/10/53	spring unit, white, strong, max. 10° range of motion				
9	SF5807-C/10/80	spring unit, yellow, very strong, max. 10° range of motion				
9	SF5807-C/05/99	spring unit, red, extra strong, max. 5° range of motion				
9	SF5807-C/0	static dorsiflexion stop, black, 0° range of motion				
System Stirrups, System Anchors and System Side Bars						
10	BR1413-L070	sliding bushing	wear	replacing	every 6 months	18 months
11	see catalogue	system stirrup	wear or breakage	replacing	every 6 months	48 months
12	see catalogue	system anchor and system side bar	wear or breakage	replacing	every 6 months	if required
-	-	adhesive bond (system side bar/system anchor and system joint)	breakage	adhering new carbon system component	every 6 months	if required

* depending on the assessment of the distributor of the custom-made product regarding the patient's usage behaviour

WP0130-GB-2026-09