
patient commission/configuration ID _____

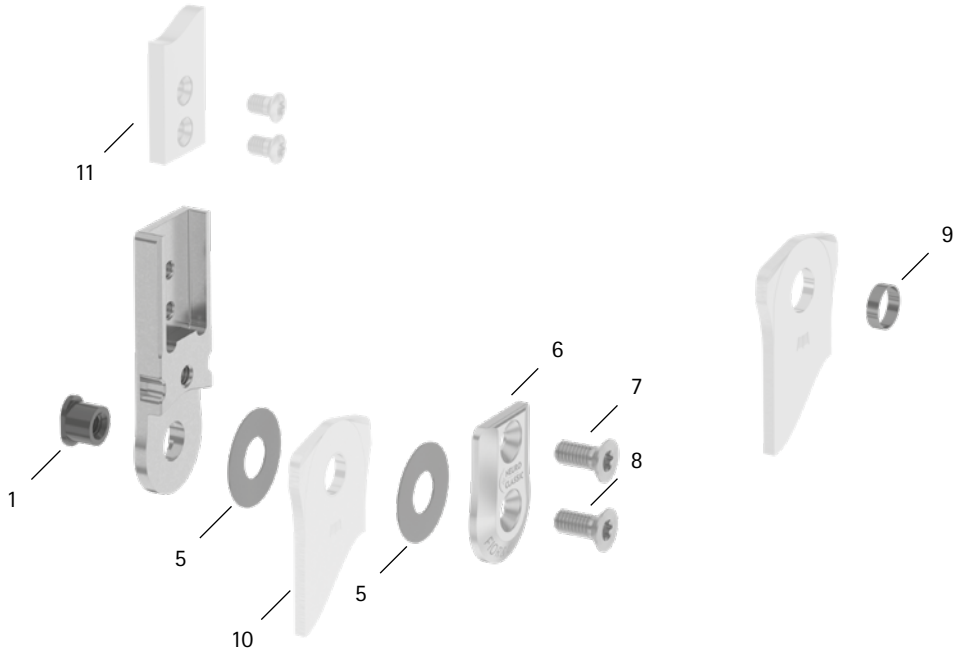

Item	24mm	Description	Potential Problem	Measure	Recommended Inspection, Potential Replacement*	Latest Replacement after
System Joint						
1	SB1262-L1190	bearing nut, short**	wear	replacing	every 6 months	36 months
5	GS3013-040	sliding washer, Ø = 30mm	wear	replacing	every 6 months	18 months
5	GS3013-045					
5	GS3013-050					
5	GS3013-055					
5	GS3013-060					
6	SF0167-3/AL	cover plate**	wear	replacing	every 6 months	36 months
7	SC1416-L16	countersunk flat head screw with hexalobular socket**	wear	replacing	every 6 months	36 months
8	SC1416-L16	countersunk flat head screw with hexalobular socket (axle screw)**				
1, 6-8	SF1987-AL	functional unit plug + go modularity***	loss of function, wear	replacing	every 6 months	36 months
System Stirrups, System Anchors and System Side Bars						
9	BR1413-L070	sliding bushing	wear	replacing	every 6 months	18 months
10	see catalogue	system stirrup	wear or breakage	replacing	every 6 months	48 months
11	see catalogue	system anchor and system side bar	wear or breakage	replacing	every 6 months	if required

* depending on the assessment of the distributor of the custom-made product regarding the patient's usage behaviour
** is part of the functional unit
*** included system components can be exchanged separately

WP0134-GB-2026-09