




Towelling Padding Material

... for comfortably padding orthoses



Towelling Padding Material


... for comfortably padding orthoses

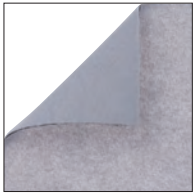

The soft padding material has a skin-friendly towelling coating on the one side and a velour coating on the reverse. The padding can be x-changed and is available in black, blue, beige and grey.



Because it can be washed at 60 °C and has an antimicrobial coating, this padding enables the best possible wearing comfort. It is a breathable, fast-drying material.

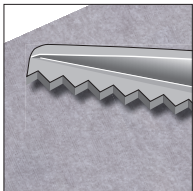

When you cut the material with pinking shears there are no frayed edges, that means you get a functional edge which does not need to be linked.



This fine material distinguishes by soft padding properties, a good resilience as well as a good adjustment to contours and offers a high wearing comfort.



To fix the padding material at its velour side, the 1mm thick, self-adhesive hook tape can be ordered separately.



A secure fixation of the padding material is achieved by selectively attaching the hook tape to the orthosis.



For hygienic reasons, we recommend to pre-cut one x-change padding at least. Thus, the patient always has a clean padding at hand.

Towelling Padding Material				
Article number	Scope of delivery	Length x width x thickness [mm]	Colour	Unit
PL3697-02/1	1 x	1000 x 1170 x 4	Black	Sheet
PL3697-02/2	1 x	1000 x 1170 x 4	Blue	Sheet
PL3697-02/3	1 x	1000 x 1170 x 4	Beige	Sheet
PL3697-02/4	1 x	1000 x 1170 x 4	Grey	Sheet



Hook tape, self-adhesive, black			
Article number	Scope of Delivery	Length x width x thickness [mm]	Unit
KV4050-L1000	1 x	1000 x 50 x 1	Band

! We recommend other materials for padding knee orthoses. We would like to give you some advice and support on this matter as well.