




# Towelling Padding Material

... for a comfortable padding of orthoses



# Towelling Padding Material


... for a comfortable padding of orthoses

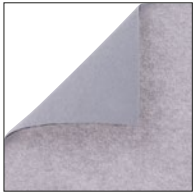

The soft padding material has a skin-friendly towelling coating on one side and a velour coating on the reverse side. It is exchangeable and available in the colours black, blue, beige and grey.



The antimicrobial properties and washability at 60°C ensure an optimal wearing comfort. The material is breathable and dries quickly.

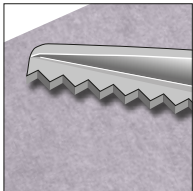

If you cut the material with pinking shears, there are no frayed edges and therefore a functional edge which does not need to be linked.



With its low thickness, the material is defined by a soft padding, good resilience and contour adaptation, thus offering a high wearing comfort.



For fixing the padding material on the velour side, hook dots or hook tape can be ordered separately.



A secure fixation of the padding material is achieved by a point fixation of the hook dots or the hook tape in the orthosis.



For hygienic reasons, we recommend producing at least one replacement padding. This way, the patient always has a clean padding at hand.

## Towelling Padding Material

Article Number	Scope of Delivery	Length x Width x Thickness [mm]	Colour	Unit
PL3687-02/1	1 x	1000 x 1400 x 4	black	sheet
PL3687-02/2	1 x	1000 x 1400 x 4	blue	sheet
PL3687-02/4	1 x	1000 x 1400 x 4	grey	sheet



## Hook Tape, Self-adhesive, Transparent

Article Number	Scope of Delivery	Length x Width x Thickness [mm]	Unit
KV4050-L1000	1 x	1000 x 50 x 1	tape

## Hook Dots

Article Number	Scope of Delivery	Unit
KV0020-22/06	6 x hook dots, self-adhesive, transparent	set
KV0020-22/40	40 x hook dots, self-adhesive, transparent	set

! For the padding of knee orthoses we recommend other materials. Please contact us for further information.